

# NEW ADVANTAGES: HIGHER TECH, LOW-PROFILE DESIGN, UNSURPASSED COMFORT.

## IMPROVED BREATHABILITY AND SKIN TOLERABILITY

**BreathSkill** is the fabric that is usually next to the skin on the supports and braces. It has been specially created to let the skin breathe, especially in warmer weather. **BreathSkill** has also been scientifically tested to ensure that it is safe for the skin and the environment.

## NEW TARGETED COMPRESSION

Most supports and braces are made with **FlexSkill**, **BreathSkill** and **CompresSkill**; three fabrics working in synergy to provide support only where needed and a gentle fit everywhere else.

## ELASTIC MEMORY

Fabrics, straps, inserts and fastenings deliver the same high level of performance over time. They lose none of their elasticity or resilience, ensuring long-lasting efficacy.

## THE BEST POSSIBLE FIT

Even easier to put on and more comfortable to wear, even for long periods of time. Even more breathable, stronger, and softer thanks to the special flocking on the **BreathSkill** fabric. Because the healing process is easier when you feel good and comfortable.



# DETAILS DESIGNED TO MAKE THE DIFFERENCE:

A patented **symmetrical fastening system** with a new **low-profile buckle** (08, 08A, 15, and 15A knee supports).



While the symmetrical fastening makes it easier to correctly position the support and prevents it from shifting when the patient moves, the new buckle made from a special polymer combines high performance with a low-profile design. Ultralight and not remotely cumbersome, it can turn horizontally and tilt vertically, perfectly adapting the circumference of any shape of leg.

## **Easyskill: the new fastening system for knee supports (open models)**

Open knee supports are fastened by overlapping and attaching the two flaps so that they do not need to be pulled on over the foot. But the larger the overlapped surface, the greater the force required to fasten the support. This is why we have redesigned and “reduced” the flaps. Smaller, more practical, more comfortable.



## Smart strap fastenings: (22, 24, 24A, 26, 26A, 28A models)



This is the new strap fastening system: just pass the straps through the holes provided, tighten, and fasten with Velcro®. Practical and smart, it retains all the adjustments and settings that were made the first time the support was fitted.

## ClickClack: the new super-fast and screwdriver-free extension and flexion adjustment system

Some knee supports feature the new extension and flexion adjustment system. Just open the hinge by sliding the two covers to the side, remove the stops holding the wedges in place and replace these wedges to provide the degree of flexion or extension prescribed by the physician. Then just secure the wedges with the stops and slide the covers back into place until they “click.”



# GENUSKILL 60 / 63

## GENUSKILL 60, SHORT TUBULAR KNEE SUPPORT

CODE REF.33160

## GENUSKILL 63, SHORT TUBULAR KNEE SUPPORT WITH OPEN KNEECAP

CODE REF.33163



- **BreathSkill**, the **3D** fabric with a woven mesh covering
- Shorter design
- (Only GenuSkill 63) Open kneecap
- **CompresSkill**, the stretchy multi-fibre yarn for compression

Code	REF.33160 / REF.33163					
Size	XS	S	M	L	XL	XXL
Thigh circ. taken 15 cm above mid-patella cm	33/39	39/45	45/52	52/58	58/66	66/75
Length cm	27	27	27	27	27	27
Colour	grey / Fits right and left					

### INDICATIONS:

- Gonalgia from bruising
- Early stage arthrosis
- (GenuSkill 60 only): 1st grade sprains to collateral ligaments
- (GenuSkill 60 only): Tendinitis
- (GenuSkill 63 only) Femoro-patellar chondropathy and chondromalacia in the knee
- (GenuSkill 63 only): Prevention of displacement of the patella during sporting activity



# genuSKILL 04

CODE REF.33104

## TUTUBLAR KNEE SUPPORT WITH PATELLA STABILISER



- **CompresSkill**, the stretchy multi-fibre yarn for compression
- **BreathSkill**, the **3D** fabric with a woven mesh covering
- **FlexSkill** fabric bonded with a technical fabric
- Patella stabiliser
- Lateral and medial stabilising stays

Code  
Size  
Thigh circ. taken 15 cm  
above mid-patella cm  
Length cm  
Colour

REF.33104

XS	S	M	L	XL	XXL
33/39	39/45	45/52	52/58	58/66	66/75
34	34	34	34	34	34

grey / Fits right and left

### INDICATIONS:

- Gonalgia from bruising
- Early stage arthrosis
- Femoro-patellar chondropathy and chondromalacia in the knee
- Mild instability of the patella
- Prevention of displacement of the patella during sporting activity



# GENUSKILL 08 / 08A

## TUBULAR KNEE SUPPORT WITH SPIRAL STAYS AND PATELLA STABILISER

CODE REF.33108

## WRAPAROUND KNEE SUPPORT WITH SPIRAL STAYS

CODE REF.33108A



- (Only GenuSkill08A) An easier-to-fasten Velcro® opening
- (Only GenuSkill08) **FlexSkill** fabric bonded with a technical fabric
- (Only GenuSkill08) Patella stabiliser
- **BreathSkill**, the **3D** fabric with a woven mesh covering
- Lateral and medial stabilising stays
- Top and bottom stabilising straps with the symmetrical fastening (rotation-proof)
  - A jointed buckle that automatically adapts to fit the leg
  - A rubber tab makes it easier to undo Velcro® fastenings
- **CompresSkill**, the stretchy multi-fibre yarn for compression

Code Size Thigh circ. taken 15 cm above mid-patella cm Length cm Colour	REF.33108 / REF.33108A					
	XS	S	M	L	XL	XXL
	33/39	39/45	45/52	52/58	58/66	66/75
	34	34	34	34	34	34
	grey / Fits right and left					

### INDICATIONS:

- Initial gonalgia after sprain trauma with mild residual collateral ligament slackness (knee inoperable)
- (GenuSkill 08 only) Mild latero-medial instability of the patella
- (GenuSkill 08 only) Femoro-patellar arthrosis



# GENUSKILL 15 / 15A

## TUBULAR KNEE SUPPORT WITH HINGED STAYS AND PATELLA STABILISER

CODE REF.33115

## WRAPAROUND KNEE SUPPORT WITH HINGED STAYS

CODE REF.33115A



- (Only GenuSkill15A) An easier-to-fasten Velcro® opening
- **BreathSkill**, the **3D** fabric with a woven mesh covering
- (Only GenuSkill 15) **FlexSkill** fabric bonded with a technical fabric
- (Only GenuSkill 15) Patella stabiliser
- Lateral and medial stabilising polycentric stays (with the option of restricting flexion)
- Top and bottom stabilising straps with the symmetrical fastening (rotation-proof)
  - A jointed buckle that automatically adapts to fit the leg
  - A rubber tab makes it easier to undo Velcro® fastenings
- **CompresSkill**, the stretchy multi-fibre yarn for compression

Code	REF.33115/ REF.33115A					
Size	XS	S	M	L	XL	XXL
Thigh circ. taken 15 cm above mid-patella cm	33/39	39/45	45/52	52/58	58/66	66/75
Length cm	34	34	34	34	34	34
Colour	grey / Fits right and left					

### INDICATIONS:

- Gonarthrosis of the after-effects of distortions with residual collateral ligament slackness (knee inoperable)
- (GenuSkill 15 only) Mild latero-medial instability of the patella
- (GenuSkill 15 only) Femoro patellar-arthritis



# GENUSKILL 22 / 22A

## TUBULAR KNEE BRACE WITH POLYCENTRIC HINGED BARS AND F-E CONTROL

CODE REF.33122

## KNEE BRACE WITH POLYCENTRIC HINGED BARS AND F-E CONTROL (OPEN)

CODE REF.33122A



- (Only GenuSkill 22A)  
An easier-to-fasten Velcro® opening
- (Only GenuSkill 22) **FlexSkill** fabric bonded with a technical fabric
- **BreathSkill**, the **3D** fabric with a woven mesh covering
- The polycentric hinge can be easily adjusted with pins on the exterior
- (Only GenuSkill 22) Patella stabiliser
- A new strap fastening system with a rapid-release buckle, ring-and-loop system, and Velcro® fastenings.
- **CompresSkill**, the stretchy multi-fibre yarn for compression

Code Size	REF.33122 / REF.33122A					
	XS	S	M	L	XL	XXL
Thigh circ. taken 15 cm above mid-patella cm	33/39	39/45	45/52	52/58	58/66	66/75
Length cm	34	34	34	34	34	34
Colour	grey / Fits right and left					

### INDICATIONS:

- Functional treatment of sprain injuries grade I and II of the collateral medial/lateral ligament
- Prevention of relapsing MCL and LCL sprains during sports activities
- Chronic latero-medial instability
- Femoro-patellar arthrosis
- Arthritic instability (inoperable knee)





# genusKILL 24A CODE REF.33124A

## LONG KNEE BRACE WITH POLYCENTRIC HINGED BARS AND F-E CONTROL (OPEN)



- **CompresSkill**, the stretchy multi-fibre yarn for compression
- **BreathSkill**, the **3D** fabric with a woven mesh covering
- An easier-to-fasten Velcro® opening
- The polycentric hinge can be easily adjusted with pins on the exterior
- A new strap fastening system with a rapid-release buckle, ring-and-loop system, and Velcro® fastenings.

**Code**  
**Size**  
**Thigh circ. taken 15 cm above mid-patella cm**  
**Length cm**  
**Colour**

REF.33124A	XS	S	M	L	XL	XXL
Thigh circ. taken 15 cm above mid-patella cm	33/39	39/45	45/52	52/58	58/66	66/75
Length cm	41	41	41	41	41	41
Colour	grey / Fits right and left					

### INDICATIONS:

- Functional treatment of sprain injuries grade I and II of the collateral medial/lateral ligament
- Prevention of relapsing MCL and LCL sprains during sports activities
- Chronic latero-medial instability
- Femoro-patellar arthrosis
- Arthritic instability (inoperable knee)



# GENUSKILL 26 / 26A

## SHORT TUBULAR KNEE BRACE WITH POLYCENTRIC HINGED BARS AND F-E CONTROL

CODE REF.33126

## SHORT KNEE BRACE WITH POLYCENTRIC HINGED BARS AND F-E CONTROL (OPEN)

CODE REF.33126A



- (Only GenuSkill 26A) An easier-to-fasten Velcro® opening
- (Only GenuSkill 26) **FlexSkill** fabric bonded with a technical fabric
- **BreathSkill**, the **3D** fabric with a woven mesh covering
- A new ClickClack polycentric hinge with wedges for easy adjustment (no screwdriver required)
- (Only GenuSkill 26) Patella stabiliser
- A new strap fastening system with a rapid-release buckle, ring-and-loop system, and Velcro® fastenings.
- **CompresSkill**, the stretchy multi-fibre yarn for compression

Codice	REF.33126 / REF.33126A					
Size	XS	S	M	L	XL	XXL
Thigh circ. taken 15 cm above mid-patella cm	33/39	39/45	45/52	52/58	58/66	66/75
Length cm	34	34	34	34	34	34
Colour	grey / Fits right and left					

### INDICATIONS:

- Functional treatment of sprain injuries grade I and II of the collateral medial/lateral ligament
- Prevention of relapsing MCL and LCL sprains during sports activities
- Chronic latero-medial instability
- Femoro-patellar arthrosis
- Arthritic instability (inoperable knee)



# genusKILL 28A

CODE REF.33128A

## LONG KNEE BRACE WITH POLYCENTRIC HINGED BARS AND F-E CONTROL (OPEN)



- **CompresSkill**, the stretchy multi-fibre yarn for compression
- **BreathSkill**, the 3D fabric with a woven mesh covering
- An easier-to-fasten Velcro® opening
- A new ClickClack polycentric hinge with wedges for easy adjustment (no screwdriver required)
- A new strap fastening system with a rapid-release buckle, ring-and-loop system, and Velcro® fastenings.

**Code**  
**Size**  
**Thigh circ. taken 15 cm above mid-patella cm**  
**Length cm**  
**Colour**

REF.33128A

XS	S	M	L	XL	XXL
33/39	39/45	45/52	52/58	58/66	66/75
41	41	41	41	41	41

grey / Fits right and left

### INDICATIONS:

- Functional treatment of sprain injuries grade I and II of the collateral medial/lateral ligament
- Prevention of relapsing MCL and LCL sprains during sports activities
- Chronic latero-medial instability
- Femoro-patellar arthrosis
- Arthritic instability (inoperable knee)



# genusSKILL 16

CODE REF.33116

## PATELLA STABILISING KNEE BRACE



- **CompresSkill**, the stretchy multi-fibre yarn for compression
- **BreathSkill**, the **3D** fabric with a woven mesh covering
- Patella alignment system:
  - lateral J-shaped foam
  - support and pressure for the kneecap with hook-and-loop straps (adjustable)
- Lateral and medial stabilising stays

<b>Code</b>	REF.33116					
<b>Size</b>	XS	S	M	L	XL	XXL
<b>Thigh circ. taken 15 cm above mid-patella cm</b>	33/39	39/45	45/52	52/58	58/66	66/75
<b>Length cm</b>	34	34	34	34	34	34
<b>Colour</b>	grey / specify right or left					

### INDICATIONS:

- Lateral-medial instability of the patella
- Femoro-patellar arthrosis
- Gonarthrosis

# malleoSkill 37

CODE REF.33137

## ANKLE SUPPORT

- **BreathSkill**, the **3D** fabric with a woven mesh covering
- **CompresSkill**, the stretchy multi-fibre yarn for compression
- Open heel

Code	REF.33137			
Size	S	M	L	XL
Ankle circumf. cm	15/20	20/25	25/30	30/36
Colour	grey / Fits right and left			

## INDICATIONS:

- 1st degree sprains of the tibiotarsal
- Tendinopathy of the instep and peronei



# malleoSkill 81

CODE REF.33181

## ANKLE SUPPORT WITH CROSSED ELASTIC STRAPS AND SPIRAL SPLINTS

- **BreathSkill**, the **3D** fabric with a woven mesh covering
- **CompresSkill**, the stretchy multi-fibre yarn for compression
- Open heel
- Two spiral support splints
- Crossover instep straps for maximum stability.

Code	REF.33181			
Size	S	M	L	XL
Circ. of ankle from base of heel across the malleoli cm	26/29	29/32	32/35	35/38
Colour	grey / Fits right and left			

## INDICATIONS:

- Conservative treatment after distortions grade II
- Acute tendinopathies
- Postoperative rehabilitation
- Supports the resumption of sporting activities after distortions grade II and III
- Chronic joint pain caused by ankle osteoarthritis



# mioSKILL 33

CODE REF.33133

## THIGH SUPPORT

- **BreathSkill**, the **3D** fabric with a woven mesh covering
- **CompresSkill**, the stretchy multi-fibre yarn for compression

Code	REF.33133				
Size	S	M	L	XL	XXL
Circumf. middle of thigh cm	48/52	52/56	56/61	61/66	66/72
Colour	grey / Fits right and left				

## INDICATIONS:

- Muscular contractions and cramps of the thigh
- Aid to resuming sporting activity after trauma
- Prevention of bruising



# mioSKILL 34

CODE REF.33134

## CALF SUPPORT

- **BreathSkill**, the **3D** fabric with a woven mesh covering
- **CompresSkill**, the stretchy multi-fibre yarn for compression

Code	REF.33134			
Size	S	M	L	XL
Circumference calf cm	33/35	35/38	38/43	43/48
Colour	grey / Fits right and left			

## INDICATIONS:

- Muscular contractions and cramps of the calf
- Aid to starting sporting activity again after trauma



## ELBOW SUPPORT

- **BreathSkill**, the **3D** fabric with a woven mesh covering
- **CompresSkill**, the stretchy multi-fibre yarn for compression
- A Velcro® strap with a nylon keeper (for increased pressure on the forearm)

Code  
Size  
Proximal circumf.  
of forearm cm  
Colour

REF.33132	XS	S	M	L	XL
	20/23	23/25	25/30	30/35	35/40
	grey / Fits right and left				

### INDICATIONS:

- Epicondylitis and epitrochleitis of the elbow
- Chronic pain from bursitis or tendinitis




# lumboSKILL 70 / 72

## LOW CORSET

CODE REF.33170

## HIGH CORSET

CODE REF.33172

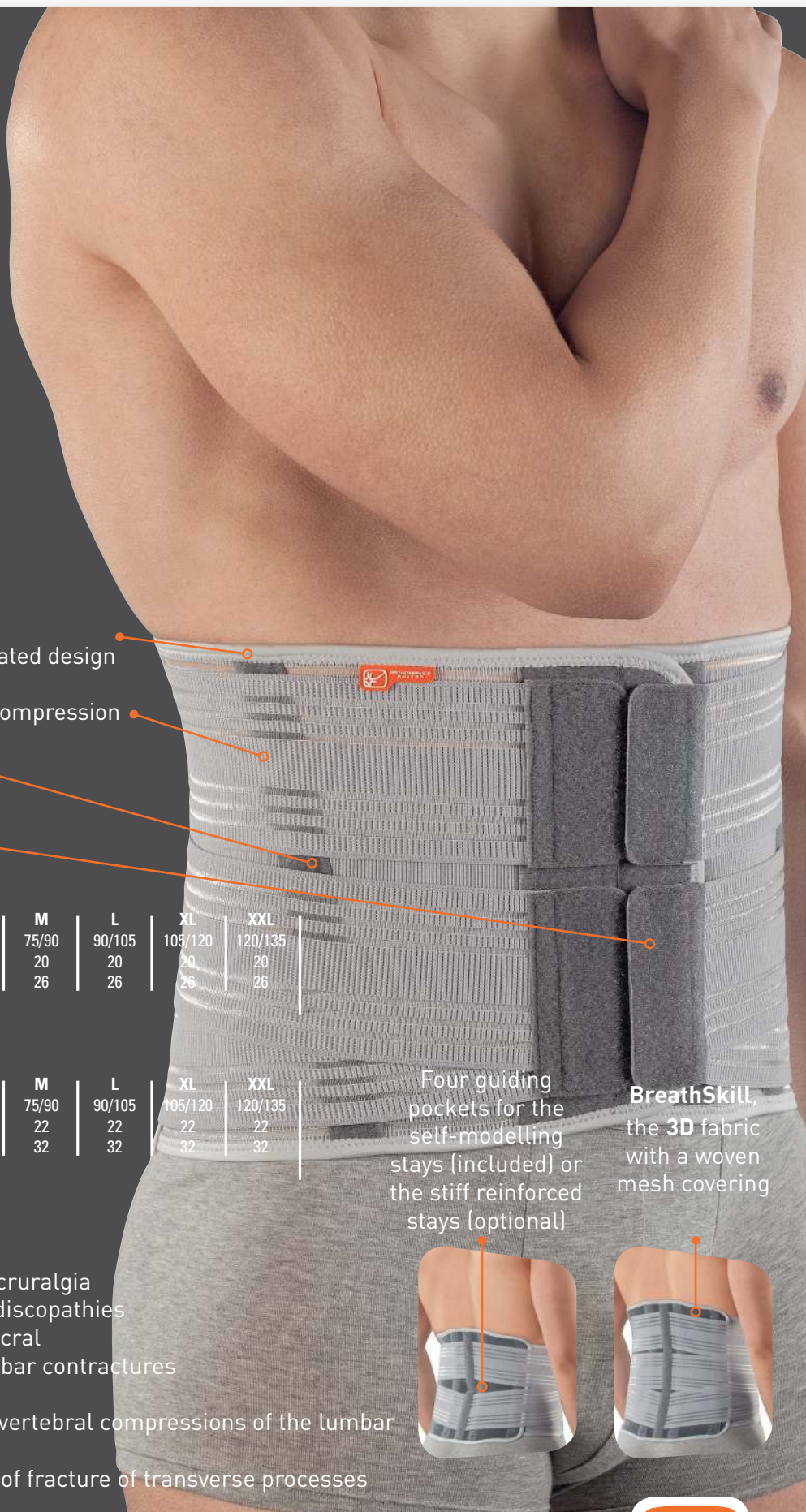

- The polyester and elastomer front and sides feature a graduated design
- High-coverage straps provide compression
- Angled spiral support stays
- Velcro® fastening

Code	REF.33170					
Size	XS	S	M	L	XL	XXL
Circumf. waist cm	50/60	60/75	75/90	90/105	105/120	120/135
Front height cm	20	20	20	20	20	20
Back height cm	26	26	26	26	26	26
Colour	grey					

Code	REF.33172					
Size	XS	S	M	L	XL	XXL
Circumf. waist cm	50/60	60/75	75/90	90/105	105/120	120/135
Front height cm	22	22	22	22	22	22
Back height cm	32	32	32	32	32	32
Colour	grey					

Four guiding pockets for the self-modelling stays (included) or the stiff reinforced stays (optional)

**BreathSkill**, the 3D fabric with a woven mesh covering

### INDICATIONS:

- Lumbalgia
- Sciatic lumbalgia and lumbar cruralgia
- Spondylarthrosis and lumbar discopathies
- Slight traumas to the lumbosacral column and paravertebral lumbar contractures
- Postoperative recovery
- LumboSkill72 only: Results of vertebral compressions of the lumbar spine
- LumboSkill72 only: Treatment of fracture of transverse processes

