

Click Medical[®] Product Catalog

Click Medical products are patent protected. For full list see www.clickmedical.co/patents

ALFAcare

Revo **Surface**[®]

Revo **Fit**[®]

Revo **Lock**[®]

Quick **Fit**[®]

Replacement Parts







Revo Surface[®]

Made with the Click[®] Reel

SKU	Product	Description	Picture
RS-100-07-01	RevoSurface [®] Tool Kit with Click [®] Reel	<p>This single RevoSurface Tool Kit with the Click Reel allows you to test adjustable designs and fit of check sockets, 3D printed devices, and on orthoses. It can also be used to solve fit issues on existing rigid devices.</p> <p>The Click Reel innovation creates adjustable devices with a lower streamlined profile, micro/macro adjustment that incrementally tightens and loosens with increased durability.</p> <p>Kit includes: 1x Click Reel, 1x Surface Collar, 1x Surface Shim, 2m HC Lace, 1.8m Lamination Tubing, 1x Metal Lace Feeder, 3x Plastic Lace Feeders, 1x Reel Tool, 10x Surface Guides, 5x Releasing Guides, 4x Phillips Head Screws, 4x Torx Head Screws</p> <p>User Weight Limit: 300lb (135kg)</p>	
RP-001	HC Lace Kit	2m of HC (High Capacity) Lace and Click Reel components for RevoSurface Tool Kit, RevoLock Align, and RevoLock 4-Hole Kits with the Click Reel.	
RP-000-07	Click [®] Reel and HC Lace Kit	<p>Replacement Click Reel and 2m of HC (High Capacity) Lace for RevoLock Align, RevoLock 4-Hole and RevoSurface Tool Kits only.</p> <p>This is not for replacing a BOA[®] dial. Please see RP0100-320-05 for replacing a BOA[®] dial.</p>	
RP-101	Surface Guide 10-Pack	The 10-Pack of Surface Guides allows a practice to have on hand additional guides for easy replacement and configuration changes. For RevoSurface Tool Kits with the Click Reel only.	



Revo Surface[®]

Made with the Click[®] Reel

SKU	Product	Description	Picture
RP-102	Surface Collar	The Click Reel mounting collar. For surface mount only. Not for use with lamination.	
RP-105	Surface Releasing Tab 5-Pack	5-Pack Surface Releasing Tabs for converting a device from a fixed lace path to a releasing lace path to fully open the device or to replace existing Surface Releasing Tabs.	
RP-201	Lamination Tube - 1.8m	Replacement Lamination Tube (1.8m) for all Revo Kits with the Click Reel and BOA [®] dials.	
BL-001	HC Lace - 25m Spool	<p>25m spool of HC (High Capacity) Lace assures you will always have replacement lace on hand.</p> <p>Cost effective, sustainable and less waste.</p> <p>Used in the RevoSurface and RevoLock Kits with the Click Reel only.</p>	
AC-100-05	Comfort Grip Ring - Black 5-Pack	Provide a soft-touch feel and 10% increased turning power to the existing Click Reel. Great for patients with dexterity issues.	
AC-101-05	Power Grip Ring - Black 5-Pack	Provide a soft-touch feel and 20% increased turning power to the existing Click Reel. Great for patients with dexterity and hand strength issues.	

Revo Fit[®]

Made with the Click[®] Reel

SKU	Product	Description	Picture
RF-200-07-01	RevoFit [®] Lamination Kit with Click [®] Reel	<p>This single RevoFit Lamination Kit with the Click Reel allows you to fabricate user-controlled adjustability anywhere on a device.</p> <p>The Click Reel innovation creates adjustable devices with a lower streamlined profile, micro/macro adjustment that incrementally tightens and loosens with increased durability.</p> <p>Kit includes: 1x Click Reel, 1x Lamination Collar, 2m HD Lace, 1x Metal Lace Feeder, 3x Plastic Lace Feeders, 1x Reel Tool, 1.8m Lamination Tubing</p> <p>User Weight Limit: 300lb (135kg)</p>	
RF-200-07-10	RevoFit [®] Lamination Fabricators 10-Pack with Click [®] Reel	<p>The RevoFit Lamination Fabricators 10-Pack with the Click Reel is recommended for central fabricators and practices that create large volumes of adjustable devices.</p> <p>Save 10% on purchase price of each single kit while assuring you will always have a kit on hand.</p>	<p>Contents:</p> <ul style="list-style-type: none"> ○ 10 Click[®] Reels ○ 10 Lamination Collars ○ 25m Spool of HD Lace ○ 5 Metal Lace Feeders ○ 20 Plastic Lace Feeders ○ 3 Reel Tools ○ 10 Lamination Tubes - 1.8m
RP-002	HD Lace Kit	2m of HD (High Durability) Lace and Click Reel components for RevoFit Lamination Kits with the Click Reel.	
RP-002-07	Click [®] Reel and HD Lace Kit	<p>Replacement Click Reel and 2m of HD (High Durability) Lace for RevoFit Kits only.</p> <p>Note: This is not for replacing a BOA[®] dial. Please see RP0100-320-05 for replacing a BOA[®] dial.</p>	

Revo Fit[®]

Made with the Click[®] Reel

SKU	Product	Description	Picture
RP-200	Lamination Collar	Replacement Lamination Collar for RevoFit & RevoLock Lamination Kits with the Click Reel only.	
RP-201	Lamination Tube - 1.8m	Replacement Lamination Tube (1.8m) for all Revo Kits with the Click Reel and BOA [®] dials.	
BL-002	HD Lace - 25m Spool	<p>25m spool of HD (High Durability) Lace assures you will always have replacement lace on hand.</p> <p>Cost effective and sustainable with less waste.</p> <p>Compatible with both the RevoFit Lamination Kit with the Click Reel and RevoFit Lamination Kit with BOA[®] dial.</p>	
AC-100-05	Comfort Grip Ring - Black 5-Pack	Provide a soft-touch feel and 10% increased turning power to the existing Click Reel. Great for patients with dexterity issues.	
AC-101-05	Power Grip Ring - Black 5-Pack	Provide a soft-touch feel and 20% increased turning power to the existing Click Reel. Great for patients with dexterity and hand strength issues.	

Revo Lock®

Made with the Click® Reel

SKU	Product	Description	Picture
RL-301-07-01	RevoLock® Align Kit with Click® Reel	<p>This single RevoLock Align Kit with the Click Reel is for transferable alignment fabrication for both upper and lower extremity sockets.</p> <p>The Click Reel innovation offers a lower streamlined profile, micro/macro adjustment that incrementally tightens and loosens with increased durability. The RevoLock Slider enhances safety for limbs changing in position or fluctuating in volume.</p> <p>Kit includes: 1x Click Reel, 1x Lamination Collar, 1x RevoLock Slider Lamination Kit, 2x RevoLock LE Quick Connect Inserts, 1x RevoLock LE Quick Connect Snap, 2m HC Lace, 1x Metal Lace Feeder, 3x Plastic Lace Feeders, 1x Reel Tool, 0.5m Lamination Tubing, 1x Align Housing + Fabrication Post</p> <p>User Weight Limit: 300lb (135kg)</p>	
RL-300-07-01	RevoLock® 4-Hole Kit with Click® Reel	<p>This single RevoLock 4-Hole Kit with the Click Reel is for standard alignment fabrication for lower extremity sockets and is easy to retrofit a standard 25x65mm shuttle lock.</p> <p>The Click Reel innovation offers a lower streamlined profile, micro/macro adjustment that incrementally tightens and loosens with increased durability. The RevoLock Slider enhances safety for limbs changing in position or fluctuating in volume.</p> <p>Kit includes: 1x Click Reel, 1x Lamination Collar, 1x RevoLock Slider Lamination Kit, 2x RevoLock LE Quick Connect Inserts, 1x RevoLock LE Quick Connect Snap, 2m HC Lace, 1x Metal Lace Feeder, 3x Plastic Lace Feeders, 1x Reel Tool, 0.5m Lamination Tubing, 1x 4-Hole Puck, 1x Tool Form Plug</p> <p>User Weight Limit: 300lb (135kg)</p>	
RL-310-00-01	RevoLock® Tool Form	<p>The 65mm x 25mm RevoLock Fabrication Tool Form is a required, reusable fabrication tool for RevoLock 4-Hole Kits.</p>	






Revo Lock®

Made with the Click® Reel

SKU	Product	Description	Picture
RP-001	HC Lace Kit	<p>2m of HC (High Capacity) Lace and Click Reel components for RevoLock Align, RevoLock 4-Hole, and RevoSurface Tool Kits with the Click Reel.</p>	
RP-000-07	Click® Reel and HC Lace Kit	<p>Replacement Click Reel and 2m of HC (High Capacity) Lace for RevoLock Align, RevoLock 4-Hole and RevoSurface Tool Kits only.</p> <p>This is not for replacing a BOA® dial. Please see RP0100-320-05 for replacing a BOA® dial.</p>	
RP-200	Lamination Collar	<p>Replacement Lamination Collar for RevoLock & RevoFit Lamination Kits with the Click Reel only.</p>	
RP-201	Lamination Tube - 1.8m	<p>Replacement Lamination Tube (1.8m) for all Revo Kits with the Click Reel and BOA® dials.</p>	
RP-305	RevoLock® Tool Form Plug 5-Pack	<p>The RevoLock Tool Form Plug 5-Pack is to ensure fabricators have extra plugs on hand.</p>	
RP-300	RevoLock® Slider Part Kit	<p>Replacement Slider and Lock Wire for all RevoLock Kits.</p>	




Revo Lock[®]

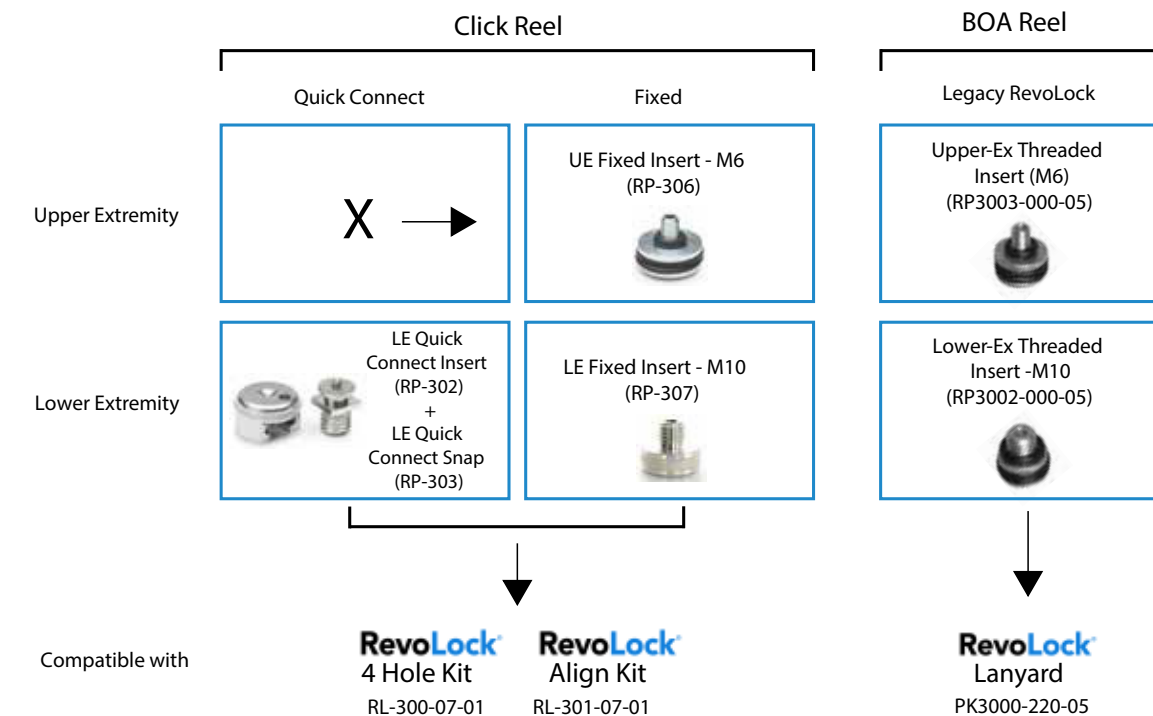
Made with the Click[®] Reel

SKU	Product	Description	Picture
RP-309	RevoLock [®] Slider Lamination Kit	Lamination Kit for the RevoLock Slider when fabricating new devices with the Click Reel.	
BL-001	HC Lace - 25m Spool	25m spool of HC (High Capacity) Lace assures you will always have replacement lace on hand. Cost effective, sustainable and less waste. Used in the RevoSurface and RevoLock Kits with the Click Reel only.	
AC-100-05	Comfort Grip Ring - Black 5-Pack	Provide a soft-touch feel and 10% increased turning power to the existing Click Reel. Great for patients with dexterity issues.	
AC-101-05	Power Grip Ring - Black 5-Pack	Provide a soft-touch feel and 20% increased turning power to the existing Click Reel. Great for patients with dexterity and hand strength issues.	
RP-306	RevoLock [®] UE Fixed Insert	RevoLock Fixed Threaded Insert (M6) for upper extremity and pediatric liners. Compatible with RevoLock 4-Hole and RevoLock Align Kits. Best for patients who cannot operate the RevoLock Quick Connect Snap or for suction systems needing an airtight seal. This part is NOT releasable and does not connect to the RevoLock LE Quick Connect Snap.	

Revo Lock[®]




Made with the Click[®] Reel

SKU	Product	Description	Picture
RP-307	RevoLock [®] LE Fixed Insert	RevoLock Fixed Threaded Insert (M10) for lower extremity liners only. Compatible with RevoLock 4-Hole and RevoLock Align Kits. Best for patients who cannot operate the RevoLock Quick Connect Snap or suction systems needing an airtight seal. This part is NOT releasable and does not connect to the RevoLock LE Quick Connect Snap.	
RP-303	RevoLock [®] LE Quick Connect Snap	Replacement RevoLock LE Quick Connect Snap connects the HC (High Capacity) Lace to the liner via the RevoLock LE Quick Connect Insert (RP-302).	
RP-302	RevoLock [®] LE Quick Connect Insert	Replacement RevoLock LE Quick Connect Insert for connection to lower extremity liners. Compatible with RevoLock LE Quick Connect Snap (RP-303) found in the RevoLock Align and RevoLock 4-Hole Kits. Often purchased for outfitting additional patient liners.	







Quick Fit[®]





Made with the BOA[®] Dial

SKU	Product	Description	Picture
OK0180-130-05	QuickFit [®] Lacer, XS	<p>The extra small size QuickFit Lacer is designed to close a small region of an orthosis, a pediatric sized orthosis, or pediatric footwear up to size 35 EU.</p> <p>Comes pre assembled with a low-profile dial and 4 guides. QuickFit Rivets sold separately, but are required to install the QuickFit Lacer kits. Lacers are sold individually. For a pair of shoes, two kits are required.</p>	
OK0185-130-05	QuickFit [®] Lacer, XS 10-Pack	<p>The extra small size QuickFit Lacer is designed to close a small region of an orthosis, a pediatric sized orthosis, or pediatric footwear up to size 35 EU.</p> <p>Comes pre assembled with 10 extra small QuickFit Lacer kits and 100 QuickFit Rivets. For a pair of shoes, two kits are required.</p>	N/A
OK0200-130-05	QuickFit [®] Lacer, SM	<p>The small size QuickFit Lacer is designed to close a small region of an orthosis, a small sized orthosis, or footwear size 36 - 41 EU.</p> <p>Comes pre assembled with a low profile dial and 6 guides. QuickFit Rivets sold separately, but are required to install the QuickFit Lacer kits. Lacers are sold individually. For a pair of shoes, two kits are required.</p>	
OK0205-130-05	QuickFit [®] Lacer, SM 10-Pack	<p>The small size QuickFit Lacer is designed to close a small region of an orthosis, a small sized orthosis, or footwear size 36 - 41 EU.</p> <p>It comes with 10 Small QuickFit Lacer kits and 100 QuickFit Rivets. For a pair of shoes, two kits are required.</p>	N/A
OK0220-130-05	QuickFit [®] Lacer, MD	<p>The medium size QuickFit Lacer is designed to close a medium sized orthosis or footwear size 42 EU and larger.</p> <p>Comes pre assembled with a low profile dial and 8 guides. QuickFit Rivets sold separately, but are required to install the QuickFit Lacer kits. Lacers are sold individually. For a pair of shoes, two kits are required.</p>	

Quick Fit[®]

Made with the BOA[®] Dial




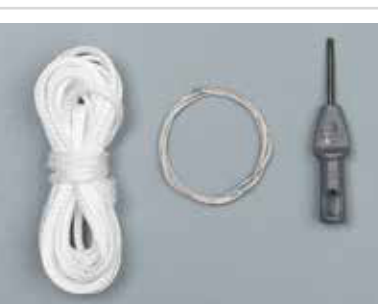
SKU	Product	Description	Picture
OK0225-130-05	QuickFit [®] Lacer, MD 10-Pack	<p>The medium size QuickFit Lacer is designed to close a medium sized orthosis or footwear size 42 EU and larger.</p> <p>It comes with 10 Medium QuickFit Lacer kits and 100 QuickFit Rivets. For a pair of shoes, two kits are required.</p>	N/A
OK0240-230-05	QuickFit [®] Lacer, LG	<p>Designed to close or secure above the ankle shoes or orthoses, such as a Leather AFO's, plastic rigid AFO's, or Cro-Mo boots. Can also be used on footwear such as hiking boots.</p> <p>Comes pre assembled with a more powerful dial and 10 guides. Rivets sold separately.</p>	
OK0165-230	QuickFit [®] Releasing Lacer, LG	<p>Adds opening functionality to our QuickFit Lacer by adding releasing guides to one side so that the system can fully open like a door for easy donning.</p> <p>Comes with a more powerful reel to close larger devices. This kit is unassembled with 4 fixed guides and 3 releasing guides. Rivets sold separately.</p>	
TL-001	QuickFit [®] Tool	<p>The QuickFit Tool is customized for use with QuickFit Rivets and is the recommended installation tool. The QuickFit Tool properly secures rivets and provides either an 1/8" and 3/16" rivet hole.</p>	
SP1040	QuickFit [®] Rivets 100-Pack	<p>QuickFit Rivets come in packs of 100. You can supply your own rivets, or purchase our easy to use rivet pack to install QuickFit Lacers to any device or shoe.</p>	

SKU	Product	Description	Picture
RP0200-130-05	QuickFit [®] XS/SM/MD Dial and Lace Replacement Kit	Replacement BOA [®] dial and 1.5m of Low-Mid Power Steel Lace for XS/SM/MD size QuickFit Lacers using BOA [®] dials.	
RP0200-230-05	QuickFit [®] LG Dial and Lace Replacement Kit	Replacement BOA [®] dial and 2m of compatible Low-Mid Power Steel Lace for LG size QuickFit Lacers using BOA [®] dials.	
LOK019	Low-Mid Power Steel Lace - 10m Spool	Bulk replacement Low-Mid Power Steel Lace (10m spool) compatible only with BOA [®] dials and QuickFit Lacer kits. NOT COMPATIBLE with Click[®] Reel.	N/A
OK1410-130-05	QuickFit [®] Strap 1.5"	<p>QuickFit Straps are designed to quickly install on a wide variety of cuffs, orthoses, and prosthetic devices to create an improved, easy to use, adjustable strap.</p> <ul style="list-style-type: none"> • Prefabricated, durable, and installs in minutes • Don/doff via simple Velcro[®] closure • 100mm (4") of adjustable travel • Usable strap length: Max: 0.75m (30") / Min: 0.2m (9") • Cut to size • Width (38mm) 	
OK1415-000-05	QuickFit [®] Buckle 1.5"	<p>The QuickFit Buckle can be added onto the QuickFit Strap in a few seconds. The magnetic buckle allows the user to easily connect the strap - even with limited dexterity or in hard to reach areas. Its plastic on plastic interface eliminates wear issues of Velcro[®] closure.</p> <ul style="list-style-type: none"> • Easily adds onto QuickFit Strap in seconds • Adds 100mm (4") to strap usable length • 5mm (.2") profile • Width (38mm) 	

Notes:



Replacement Parts

using the BOA® dial

SKU	Product	Description	Picture
PK0060-000-10	RevoFit® Re-Lacing Kit for BOA® dial	2m of HD (High Durability) Lace for RevoFit Lamination Kits Compatible only with RevoFit Kits with BOA® dial.	
RP0100-320-05	RevoFit® Dial and Lace Replacement Kit with BOA® dial	Replacement BOA® dial and 2m of HD (High Durability) Lace Compatible only with RevoFit Kits with BOA® dial.	
RP3000-220-05	RevoLock® Lanyard Dial and Lace Replacement Kit with BOA® dial	Replacement BOA® dial and 1.2m of lace Compatible only with RevoLock Kits with BOA® dial.	
RP3001-000-05	RevoLock® Lanyard Re-Lacing Kit for BOA® dial	Contains 1.2m of lace Compatible only with RevoLock Kits with BOA® dial.	

Replacement Parts

using the BOA® dial

SKU	Product	Description	Picture
RP3002-000-05	RevoLock® Lanyard - Lower-Ex Threaded Insert (10mm)	Replacement RevoLock Lower-Extremity Threaded Insert which attaches the lace to the liner via a threaded connection. NOT COMPATIBLE with LE Quick Connect Snap and insert found in the RevoLock 4-Hole or RevoLock Align Kits with the Click Reel.	
RP3003-000-05	RevoLock® Lanyard - Upper-Ex Threaded Insert (6mm)	Replacement RevoLock Upper-Extremity Threaded Insert which attaches the lace to the liner via a threaded connection. NOT COMPATIBLE with LE Quick Connect Snap and insert found in the RevoLock 4-Hole or RevoLock Align Kits with the Click Reel.	
LOK020	RevoLock® Soft Lace - 10m Spool	Bulk replacement lace (10m spool) Compatible only with RevoLock Kits with BOA® dial.	N/A

AlfaCare AS
Heddalsveien 11, bygg 140
3674 Notodden
E-post: post@alfacare.no
Telefon: 35 02 95 95